HOMEDICS[®]
B E A U T Y

Remove

MICRODERMABRASION + COOLING SURFACE

3 YEAR GUARANTEE

MD-1000-EU2


Register your product today at
www.homedics.co.uk/product-registration

PRODUCT FEATURES

HoMedics Microdermabrasion brings professional salon skincare into your own home. The diamond tip works with a gentle but powerful suction action to remove dead cells from the skin surface, revealing and stimulating the fresh skin layer beneath. Perfect for restoring youthful radiance to the skin; improving skin tone and texture, smoothing the appearance of fine lines and wrinkles and reducing the appearance of age spots and sun damage. Use regularly to maintain young and fresh looking skin.

1. Diamond tip
2. Black wool filter
3. Handle
4. Charging base
5. Power button & level selector
6. Cooling button
7. Cooling surface
8. Replacement black wool filters
9. Cleaning brush
10. Cloth travel bag
11. USB cable

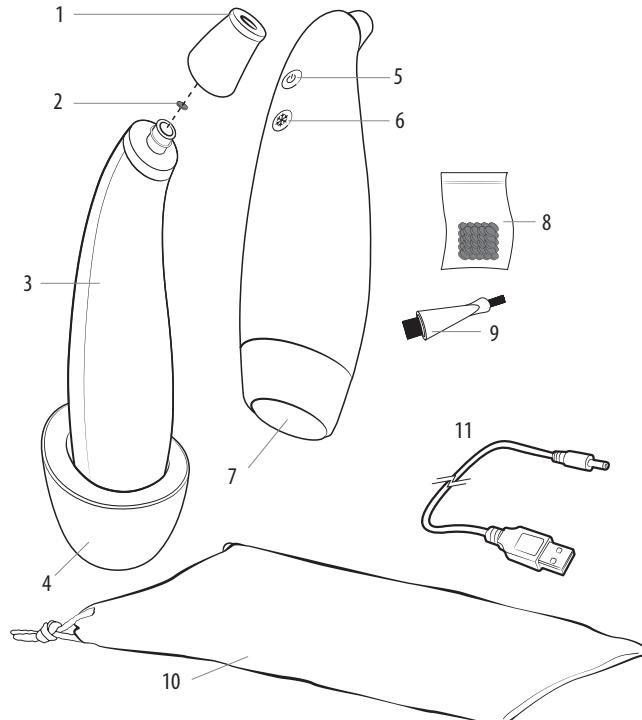UNIQUE FEATURES & SPECIFICATIONS

REMOVE Microdermabrasion + Cooling Surface brings professional spa skincare into your own home. The diamond tip works with a gentle but powerful suction action to remove dead cells from the skin's surface, revealing and stimulating the fresh skin beneath. Perfect for restoring youthful radiance to the skin, improving skin tone and texture, smoothing the appearance of fine lines and wrinkles, and reducing the appearance of age spots and sun damage. Use regularly to maintain youthful-looking skin.

CLINICALLY TESTED FOR ALL SKIN TYPES

Spa-quality microdermabrasion device

4 INTENSITY LEVELS

Adjusts for maximum comfort and effectiveness

2 DIAMOND TIPS

Normal and sensitive tips

COOLING SURFACE

Soothes skin post-treatment

BLACK WOOL FILTERS

100 black wool filters

RECHARGEABLE

Charging base and charging cord

INSTRUCTIONS FOR USE

AREA	Cheeks	Around Lips & Chin	Around Eyes	Forehead	Nose
MAX LEVEL	3	2	1	4	4
ROUTINE					

NOTE: Avoid the sensitive areas immediately above and below the eye.

Charging the battery

When the level indictor LED flashes red, the product needs charging. Connect the charging base to the mains using the USB power lead and the appropriate plug. Place the handle in the charging base and press down firmly. The LED lights will flash to indicate that charging is taking place. When charging is complete all the lights will stay on. It takes 4 hours for a full charge, which will provide 90-120 minutes use, dependent on power levels used.


During use you should feel a gentle suction effect as the tip moves across the skin and some people may also notice the slight coarseness of the diamond tip.

Depending on your skin type, you may experience some reddening or increased sensitivity after treatment. This is perfectly normal and usually subsides within a few hours. You may wish to consider this and make your first treatment at the end of the day rather than just before going out.

Adverse reactions are rare, but if skin becomes unduly inflamed or irritated, stop using the device immediately.

Filter Fitting and Replacement

For correct operation it is very important that a new black wool filter is fitted before each use of the device. The filter collects dead skin cells and should be disposed of after use and then replaced. Using the device without a filter will allow debris and dead skin cells to enter the motor, shortening the life of the product and invalidating the guarantee.

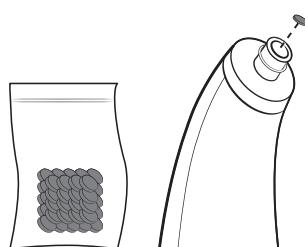What to Expect

Regular microdermabrasion treatments can both restore and maintain a youthful appearance to the skin. We advise starting slowly with just one or two treatments a week, then building up to treatments every other day if desired.

Facial Treatment Routine

- Wash and dry your face
- Remove the handle from the charging base and fit your choice of diamond tip with a new black wool filter.
- Press and hold the power button for 2 seconds to turn the unit on. Press again as required to increase the power level, as indicated by the LEDs. It is recommended to start on level 1 for the first treatment and to progress to a higher level if comfortable. See also MAX recommended levels in the above diagrams
- Place the diamond tip flat against the skin until you feel suction then move it slowly across the surface of the skin in a smooth motion, covering the same area no more than 2-3 times in one treatment. Do not press hard. In some areas it can help to hold the skin taught with the other hand (eg. around the eyes).
- The diagrams above will help you establish an effective facial routine.

After treatment cooling surface

- Press and hold the cooling button for 1 second. A blue LED will glow indicating the cooling surface is activated
- Apply the cooling surface to the treated area to soothe the skin post-treatment.

- DO NOT hold the cooling surface on the same area of skin for more than 10 seconds. The product is designed to be moved continuously across the surface of the skin.
- Press and hold the cooling  button for 1 second to turn the cooling surface off

NOTE: the cooling surface is equipped with a 2.5 minute auto-off feature.

NOTE: the cooling surface operates independently. It is normal for the unit to feel warm when using the cooling surface

After treatment

- Press and hold the power button for 2-3 seconds to turn off device.
- Wash your face with cold water to remove any remaining dead skin cells and then apply your preferred moisturiser.
- NOTE:** Avoid using AHA (acid based) moisturisers on the day of treatment
- Depending on your skin type, you may experience some reddening or increased sensitivity after treatment. This is perfectly normal and usually subsides within a few hours.
- Avoid direct sun exposure after treatment, consider applying a stronger sunscreen as required.
- Always, dispose of the used black wool filter and replace with a new filter ready for next use.
- Remove dead skin cells and debris from the diamond tip using the brush provided. We also recommend the periodic use of an ultrasonic cleaner for deep cleaning the diamond tip, especially if it starts to feel smooth or less effective.

Other Treatments

Other areas of the body can also be treated 2-3 times a week.

For example, neck, chest, hands etc. When treating larger areas, replace the black wool filter part way through the treatment in order to maintain the level of suction. Aim to cover the same area no more than 2-3 times in one treatment.

FAQ

How often should I use the product?

It's best to start slowly – usually just once a week, then build up to 2 or 3 times a week if this is comfortable for your skin.

Will it hurt?

No. Treatment should not be painful. Depending on the sensitivity of your skin, you may feel slight discomfort from the diamond tip as it moves across the skin and there could be some short term redness or sensitivity after treatment, but this will soon subside. If treatment is uncomfortable, then you may need to use a lower power setting or switch to the sensitive diamond tip.

Can I use the product in combination with anti-aging creams/serums?

If your anti-aging cream or serum contains AHA (or acid) then it should not be used on the same day as the microdermabrasion

treatment, but can be used on non-treatment days. Non-AHA (acid free) anti-aging products can be used at any time.

The product does not switch on

The power button has to be pressed continuously for 2-3 seconds to turn the product on. If you try this and it still doesn't work, then it may need recharging. See 'CHARGING THE BATTERY'

The product seems to be losing power

This can happen if the black wool filter is not replaced after each treatment. Try fitting a new filter. Also check the tip to ensure that the plastic part of the tip has no cracks or damage. Check that the two rubber o-rings are in place on the nozzle end of the unit (where the tip fits over).

The diamond tip no longer seems rough

The diamond tip can accumulate skin cells over time which may start to make it feel smooth and reduce its effectiveness. If this happens we recommend cleaning the diamond tip in an Ultrasonic Cleaner. These are either available to purchase, or jewellery stores offer a cleaning service. If using the product regularly you may need to do this several times a year to maintain optimum effectiveness.

What if I still have a problem or question that is not answered here?

Please contact our service team who will be very happy to help you.

Customer Service

Tel: +44 (0) 1732 378 557

Email: support@homedics.co.uk

www.homedics.co.uk

MAINTENANCE & STORAGE

To Store

Place attachments into their proper compartments inside the storage case. Store case in a safe, dry place, out of children's reach.

To Clean

Always unplug before cleaning the unit. Never use liquids or abrasive cleaners to clean the housing of the unit or the case. Wipe off case with a damp cloth. Clean attachments periodically using detergent. Allow attachments to dry before inserting them into the hand unit.

ACCESSORIES AND SPARE PARTS

Available from our website: www.homedics.co.uk

- Replacement Black wool filters (100pcs)
- Diamond tip
- Sensitive diamond tip
- Precision tip - for precise treatment of small areas
- Body tip - for use over larger areas.

**READ ALL INSTRUCTIONS BEFORE USE.
SAVE THESE INSTRUCTIONS FOR FUTURE
REFERENCE.**

- This appliance can be used by children aged from 14 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- DO NOT reach for an appliance that has fallen in water or other liquids. Keep dry – DO NOT operate in wet conditions.
- NEVER insert pins, metallic fasteners or objects into the appliance or any opening.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. All servicing of this appliance must be performed at an authorised HoMedics Service Centre.
- Please ensure that all hair, clothing and jewellery are kept clear of the product at all times.
- If you have any concerns regarding your health, consult a doctor before using this appliance.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP.
- Pregnant women, diabetics and individuals with pacemakers should consult a doctor before using this appliance. Not recommended for use by individuals with sensory deficiencies including diabetic neuropathy.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- Do not use on irritated, sunburned, or chapped skin, or on any area with dermatological problems
- Do not apply over warts, moles, or varicose veins.
- Do not use for longer than the recommended time.
- This product contains a rechargeable battery and should not be exposed to excessive heat. Do not leave in direct sun or near a heat source such as a fire. The battery must not be replaced by the user.
- Failure to follow the above may result in the risk of fire or injury.

- **Microdermabrasion is NOT recommended if you have any of the following conditions:**

- Lesions, warts, or recent herpes outbreak
- Sunburn or sun sores
- Active rosacea
- Auto-immune disorders
- Skin cancer
- Vascular lesions
- Any open wounds, sores, or cuts
- Taking oral blood thinners (anti coagulants)

3 YEAR GUARANTEE

FKA Brands Ltd guarantees this product from defect in material and workmanship for a period of 3 years from the date of purchase, except as noted below. This FKA Brands Ltd product guarantee does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorised accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of FKA Brands Ltd. This guarantee is effective only if the product is purchased and operated in the UK / EU. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and / or authorised, or repair of products damaged by these modifications is not covered under this guarantee. FKA Brands Ltd shall not be responsible for any type of incidental, consequential or special damages.

To obtain guarantee service on your product, return the product post-paid to your local service centre along with your dated sales receipt (as proof of purchase). Upon receipt, FKA Brands Ltd will repair or replace, as appropriate, your product and return it to you, post-paid. Guarantee is solely through HoMedics Service Centre. Service of this product by anyone other than the HoMedics Service Centre voids the guarantee. This guarantee does not affect your statutory rights.

For your local HoMedics Service Centre, go to
www.homedics.co.uk/servicecentres

Battery replacement

Your product includes a rechargeable battery designed to last the lifetime of the product. In the unlikely event that you should require a replacement battery, please contact Customer Services, who will supply details of the warranty and out-of-warranty replacement battery service.

Battery directive

 This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK
FKA Brands Ltd, 29 Earlsfort Terrace, Dublin 2, Ireland